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
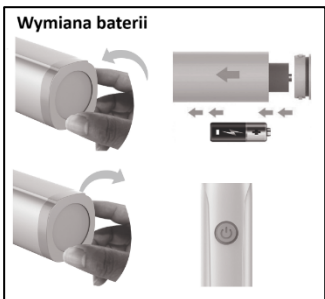
USER MANUAL U267

Importer: VITOLOG.PL Sp. z o.o, ul. Międzyńska 6, 50-514 Wrocław, Poland

Product features: Sonic toothbrush for children. Number of tips: 4. Number of operating modes: 5. Charging station: none. Waterproof: IPX7. Power supply: USB cable (500 mAh battery). Rated power: 1 W. Rated voltage: 1.5V. Input voltage: 5V. Rated voltage: 3.7 V. Color: navy blue. Material: ABS.

How to use: To turn on the device, press the switch. After turning on the device, you can change the brush's operating mode using the mode indicator. A long press on the switch turns it off. Using the power button you can change the mode and power of the device. The toothbrush has a memory function and remembers the last operating mode. The device works for 2 minutes with a break every 30 seconds to remind you to change the brushing area. When the battery level is low, the device's light flashes and the toothbrush works with less power.

Battery replacement: turn the battery compartment cover counterclockwise. Insert the battery, remembering to

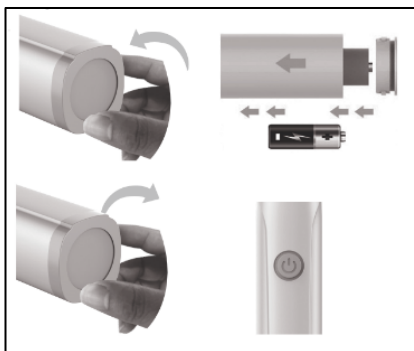observe the correct polarity. Screw the battery compartment clockwise.

Brushing your teeth: Squeeze a small amount of toothpaste onto the brush heads (less than on a regular toothbrush). Gently place the toothbrush on your toothbrush and turn it on. Instead of using the same force as a manual toothbrush, gently place the toothbrush bristles on the gum line and use high-frequency vibrations to clean. After 30 seconds, move to the next brushing area.


Replacing the head: replace the brush head every 1 month. After some time of use, the toothbrush will become rough and will easily collect bacteria, making it impossible to achieve cleaning results.

Cleaning the device: when washing, turn your head downwards.

After each use, wash the bristles and brush head and then dry them. Do not pull the brush bristles so that they do not deform and fall out. The toothbrush body should be cleaned with a dry cloth and dried. Do not immerse the toothbrush in water for a long time. Do not use water above 40°C or acidic cleaning agents for cleaning. The brush head needs to be replaced. It is recommended to replace the brush head every 1 month or when the bristles are bent or deformed (actually depends on individual use).

Warnings: Please read the instruction manual before use and keep it for future reference. Use the product as intended. The product is intended for home use. When using this product for the first time, due to the rapid vibration of the bristles, some people will feel a stimulating numbness, do not worry about it, it is a normal phenomenon. This feeling will gradually decrease with subsequent use to adapt to the vibration frequency of the toothbrush. This product can only be used with a special brush head. This product is intended only for cleaning teeth, gums and tongue. Do not use for other purposes. This product is for personal use only and is not intended for use by dental clinics or multi-patient institutions. Do not expose the toothbrush to direct sunlight. Store in a dry and ventilated place at a temperature between 7°C and 40°C. For indoor use. This product is not intended for use by persons with disabilities, mental illness or lack of knowledge or experience (including children), unless they are supervised by a guardian. After the first few uses of an electric toothbrush, you may notice slight bleeding of the gums, which usually subsides after a few days. If your gums continue to bleed, consult your dentist. If your mouth or gums have undergone surgery within the last two months, consult your dentist before using this toothbrush. Do not use hair care products (such as gel, oil) or hand sanitizers in contact with this product to avoid damaging its surface. Do not use brush heads that are not recommended by the manufacturer. Store and use the product in a dry, dust-free place, away from sources of moisture. Keep out of reach of children. Keep the product clean. Do not use caustic agents, bleaches, solvents or rough surfaces to clean. Protect the product against strong vibrations, falls and mechanical

damage. Regularly check the technical condition of the device and its charger. Do not use a damaged product. Any repairs may only be performed by an authorized person. In case of product failure, do not touch the device and disconnect it from the power source. Do not open the device. Do not modify the device or perform any repairs. If any problems occur with your device, stop using it. It is normal for the toothbrush to become deformed during production and transport. It may return to normal after soaking in warm water. To better maintain the battery, please charge it once every 10-15 days. Do not use the charger outdoors or near hot surfaces. Unplug the charger after each charging to avoid problems such as scrapping the charger or electric shock due to overheating. Do not use the charger in a humid environment, keep the charger dry when charging. Do not connect the charger with wet hands. Before use, charge the device and disconnect it from the power source. When charging the device, make sure that the input current and input voltage match the device parameters. Do not charge the device if you notice that the cable or plug is damaged. Before each subsequent use, make sure the device is properly charged. The appliance must not be used by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction. The product is not a toy. Do not allow children to play with.



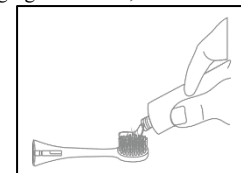
Information for users on the disposal of electrical and electronic equipment

(for households) The symbol shown on the products or the accompanying documentation indicates that faulty electrical or electronic equipment cannot be disposed of together with household waste. The correct procedure, if necessary, to dispose of electrical or electronic equipment, to dispose of, reuse or recover components, is to transfer

the equipment to a specialized collection point, where it will be accepted free of charge. In some countries, you can give the product back to your local distributor when you purchase another device. Proper disposal of the equipment preserves valuable resources and avoids negative health and environmental impacts that may be compromised by inadequate waste management. Please contact your local authority for details on the nearest collection point. Improper disposal of waste is punishable by penalties provided for in the relevant local legislation.

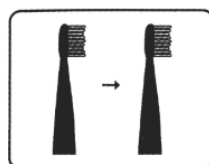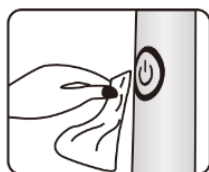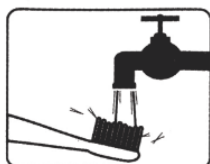Users in European Union countries: If you need to dispose of electrical or electronic equipment, please contact your nearest point of sale or your supplier for further information. **Disposal of waste in countries outside the European Union:** This symbol is only valid in the European Union. If you need to dispose of this product, please contact your local authorities or dealer for information on how to proceed correctly.

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